

SteriPure-P Polyethersulfone Pleated Cartridge Filter

»»» Sterilize grade PES filter with Longer Service Time and High Flow Rate



Specifications

» Materials of Construction

Filter membrane	Polyethersulfone
End-Filter membrane	Polyethersulfone
Supports	Polypropylene
Cage/End Caps	Polypropylene/SS 316L
Core	Polypropylene
Adaptor	Polypropylene
O-rings	Please refer to the ordering info.

» Operating Parameters

Maximum Operating Temperature	80°C
Maximum Differential Pressure (forward)	5.2bar@25°C
Maximum Differential Pressure (reverse)	2.1bar@25°C

» Filter Dimensions

Outer Diameter	69.0mm
Inner Diameter	33.0mm
Filtration Area	0.55m ² /10inch

» Safety

Bacterial Endotoxin	≤0.25EU/mL
Extractable	≤40mg/10inches

» Sterilizable

Autoclave	124°C, 30min, 30cycles
	134°C, 20min, 30cycles
Steam In-Place	124°C, 30min, 30cycles
	134°C, 20min, 30cycles

Description

SteriPure-P cartridge filters have a built-in pre-filter membrane in front of the end-filter membrane. The unique anisotropic Polyethersulfone pre-filter membrane reduces the plugging rate of end-filter membrane effectively. The optimized SteriPure-P cartridge filters can provide longer service time than generic cartridge filters in most applications.

» Reliability

From raw materials purchase, transport to production and storage, all operations follow ISO 9001 quality management system. SteriPure-P manufactured, tested and packaged in a clean room to ensure product cleanliness. SteriPure-P cartridge filters are 100% integrity tested before they shipment. A validation guide for compliance with regulatory requirements is available.

» Compatibility

SteriPure-P cartridge filters are sealed using thermal bonding process, contain no adhesive and surfactant. The components of SteriPure-P, include of Polyethersulfone, Polypropylene and Silicone/EPDM, provide broad chemical compatibility from pH 1 to pH14 and low extractable levels at high temperature.

» Economy


SteriPure-P cartridge filters take use of the asymmetric, hydrophilic high porosity polyethersulfone membrane as pre-filter and end-filter. SteriPure-P cartridge filters are specifically designed to deal with broad variety of contaminants in various applications. SteriPure-P cartridge filters can provide the highest process efficiency and minimized filtration costs.

» Regulatory Compliance

FDA 21 CFR 177
ISO 10993-Part 1, 5
EN 285:2015
Regulation (EC) No 1935/2004

» Longer Service Time and Excellent Cost-Effective


The built-in pre-filter membrane significantly increases the total throughput and stability of flow rate. In most time of your filtration process, SteriPure-P cartridge filters can provide high flow rate and low pressure drop due to the high asymmetric pre-filter membrane and controlled asymmetric end-filter membrane. SteriPure-P can effectively reduce the operating costs of filtration system, the cost of replacement filters and downtime in many applications.

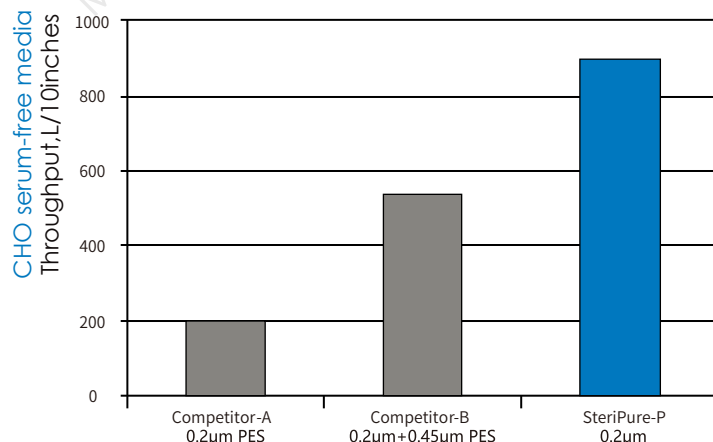

Application

- Sterilization filtration of water for injection, cleaning solutions, purified water, deionized water, etc.
- Sterilization filtration of LVP, APIs, buffers, eye drops, disinfectants, etc.
- Sterilization filtration of vaccines, serums, biological products, antibiotic water-based liquid, etc.
- Filtration of redwine, beer, juice, purified water, etc.

Soy Peptone Supplemented Culture Media
Throughput, L/10 inches

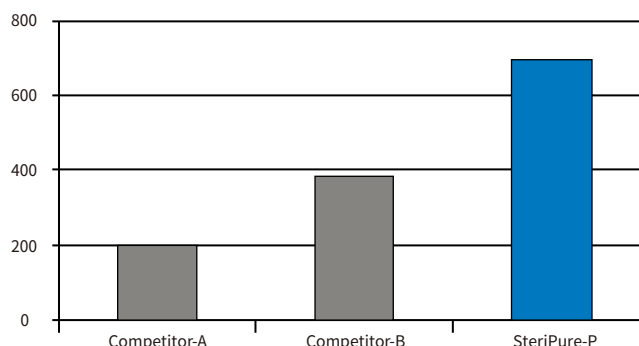

Suitable for Severe High Temperature Sterilizing


SteriPure-P can ensure integrity at temperature up to 134°C and pressure drop up to 0.5bar due to the high strength of Polyethersulfone membrane and enhanced Polypropylene structure. SteriPure-P can provide excellent filtration service for your process.

» Integrity Test Data*

Pore Size	Bubble Point	Diffusion
0.10µm+0.20µm	≥4.2bar	≤16ml/min@ 3.4bar
0.20µm+0.45µm	≥3.5bar	≤16ml/min@ 2.6bar
0.20µm+0.80µm	≥3.5bar	≤16ml/min@ 2.6bar
0.45µm+1.0µm	≥1.8bar	≤16ml/min@ 1.4bar

*Test filter cartridge length 10 inches, Pure Water @25°C

Typical Liquid Flow Rates

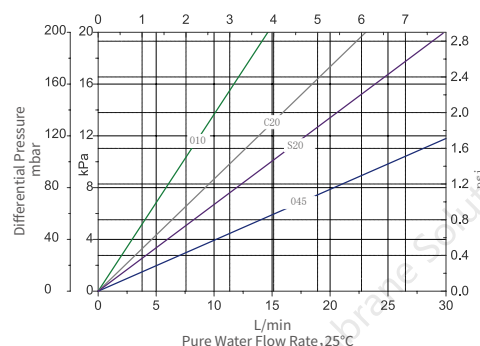

(Water / 1 cP fluid)

Typical flow rates. For liquids with viscosity differing from water, multiply the pressure drop by the viscosity in centipoise.
Unit conversion: 1 bar = 100 kPa = 14.5 psi

Ordering Information*

Filter Media	Length	Pore Size	Adaptor	Sealing	Core	Grade
CRPESP: Polyethersulfone membrane	005 = 5 inches 010 = 10 inches 020 = 20 inches 030 = 30 inches 040 = 40 inches	010=0.10µm+0.20µm 020=0.20µm+0.80µm C20=0.20µm+0.45µm (Standard) S20=0.20µm+0.45µm (High dirt holding) D45=0.45µm+0.65µm E45=0.45µm+0.80µm	1=222/Flat with SS insert 4=222/Fin with SS insert 5=226/Fin with SS insert 8=226/Flat with SS insert	S=Silicone E=EPDM N=NBR F=FKM	Blank=PP R=GF Reinforced PP S=316L Stainless Steel	PR=Pharmaceutical Requirements

* This information is a guide to the part number structure and possible options. For applications requiring autoclaving and sterilization, MS recommends the use of adaptor with SS insert to ensure filter sealing after cooling. If products that require autoclaving and sterilization use EPDM sealing rings and require customization. For availability of specific options and housing details, please contact MS.



Phone: 001-866-528-4572
Email: info@membrane-solutions.com
Address: 4530 B STREET NW #A, AUBURN WA 98001, USA



China

Phone: 021-61478109
Email: info@membrane-solutions.com
Address: 2202, No.1759 North Zhongshan Road, Shanghai 200061